


Broadband Faraday Isolator



Broadband ISO Series User's Manual

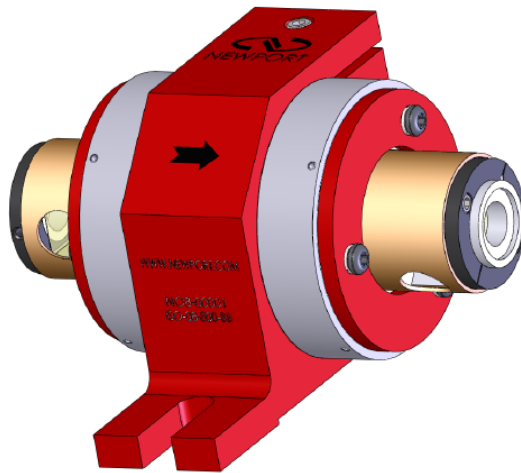

Newport®
Experience | Solutions

Warranty

Newport Corporation warrants that this product will be free from defects in material and workmanship and will comply with Newport's published specifications at the time of sale for a period of one year from date of shipment. If found to be defective during the warranty period, the product will either be repaired or replaced at Newport's option.

To exercise this warranty, write or call your local Newport office or representative, or contact Newport headquarters in Irvine, California. You will be given prompt assistance and return instructions. Send the product, freight prepaid, to the indicated service facility. Repairs will be made and the instrument returned freight prepaid. Repaired products are warranted for the remainder of the original warranty period or 90 days, whichever is longer.

Limitation of Warranty

The above warranties do not apply to products which have been repaired or modified without Newport's written approval, or products subjected to unusual physical, thermal or electrical stress, improper installation, misuse, abuse, accident or negligence in use, storage, transportation or handling. This warranty also does not apply to fuses, batteries, or damage from battery leakage.

THIS WARRANTY IS IN LIEU OF ALL OTHER WARRANTIES, EXPRESSED OR IMPLIED, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR USE. NEWPORT CORPORATION SHALL NOT BE LIABLE FOR ANY INDIRECT, SPECIAL, OR CONSEQUENTIAL DAMAGES RESULTING FROM THE PURCHASE OR USE OF ITS PRODUCTS.

First printing 2005

© 2005 by Newport Corporation, Irvine, CA. All rights reserved. No part of this manual may be reproduced or copied without the prior written approval of Newport Corporation.

This manual has been provided for information only and product specifications are subject to change without notice. Any change will be reflected in future printings.

Newport Corporation
1791 Deere Avenue
Irvine, CA, 92606 USA
P/N 90025451 Rev. A

Technical Support Contacts

North America & Asia

Newport Corporation Service Dept.

1791 Deere Ave. Irvine, CA 92606

Telephone: (949) 253-1694

Telephone: (800) 222-6440 x31694

Europe

Newport/MICRO-CONTROLE S.A.

Zone Industrielle

45340 Beaune la Rolande, France

Telephone: (33) 02 38 40 51 56

Asia

Newport Opto-Electronics Technologies

253 Aidu Road, Bld #3, Flr 3, Sec C,

Shanghai 200131, China

Telephone: +86 21 5046 2300

Telephone: +86 21 5046 2323

Newport Corporation Calling Procedure

If there are any defects in material or workmanship or a failure to meet specifications, promptly notify Newport's Returns Department by calling 1-800-222-6440 or by visiting our website at www.newport.com/returns within the warranty period to obtain a **Return Material Authorization Number (RMA#)**. Return the product to Newport Corporation, freight prepaid, clearly marked with the RMA# and we will either repair or replace it at our discretion. Newport is not responsible for damage occurring in transit and is not obligated to accept products returned without an RMA#.

E-mail: rma.service@newport.com

When calling Newport Corporation, please provide the customer care representative with the following information:

- Your Contact Information
- Serial number or original order number
- Description of problem (i.e., hardware or software)

To help our Technical Support Representatives diagnose your problem, please note the following conditions:

- Is the system used for manufacturing or research and development?
- What was the state of the system right before the problem?
- Have you seen this problem before? If so, how often?
- Can the system continue to operate with this problem? Or is the system nonoperational?
- Can you identify anything that was different before this problem occurred?

Table of Contents

Warranty	ii
Technical Support Contacts	iii
Table of Contents	iv
1 General Information	1
1.1 Introduction	1
1.2 Safe Use of Your Newport Faraday Isolator	2
1.3 Operation	4
1.3.1 The Newport Broadband Faraday Isolator	4
1.3.2 Using your Faraday Isolator	6
1.3.3 Tuning your Faraday Isolator	7
1.4 Specifications	7
2 Factory Service Information	9
2.1 Service Form	9

1 General Information

1.1 Introduction

Newport's ISO-05-800-BB and ISO-08-800-BB Faraday Isolators are essentially uni-directional light valves or spectrally selective one-way mirrors. They are used to protect a laser source from destabilizing feedback or actual damage from back-reflected light. Figure 1 below identifies the main elements of your Faraday Isolator.

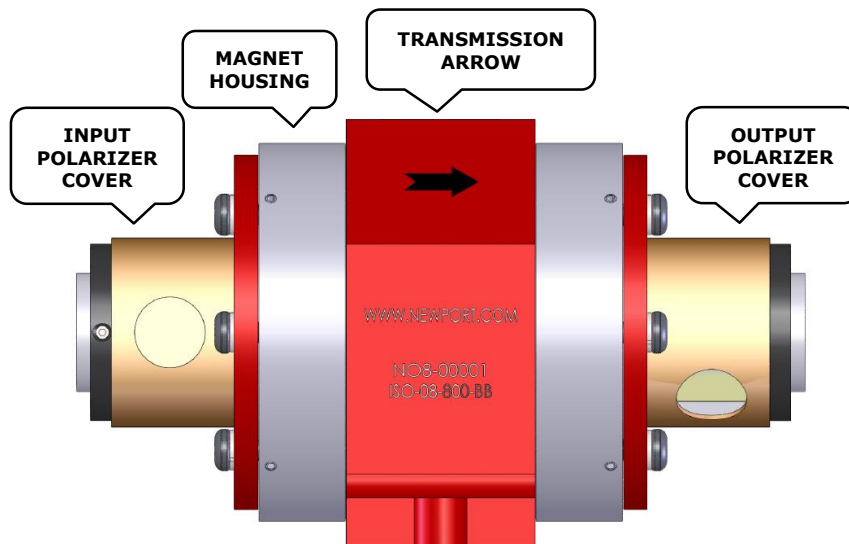


Figure 1: ISO-08-800-BB, Broadband, 8mm Aperture Faraday Isolator

The 5mm and 8mm aperture Faraday Isolators are cylindrically-shaped magneto-optic devices. Strong Neodymium Iron Boron permanent magnets are used to generate high (>10,000 Gauss) axially-oriented fields within the magnet housing. The strong longitudinal field causes 45 degrees of non-reciprocal polarization rotation for propagating light via the Faraday Effect in the Terbium Gallium Garnet ("TGG") crystal located within the magnet housing. Following the TGG crystal is a quartz reciprocal rotator with 45 degrees rotation. In the forward direction, the two rotations add up, resulting in 90 degrees of rotation. In the reverse direction, the two rotations are opposite and result in 0 degrees of rotation. The change in rotation as the wavelength shifts from the central wavelength is similar for both TGG and quartz, resulting

in broadband operation. In operation, the magnet housing is sandwiched between input and output polarizers that have their transmission axis oriented 90 degrees relative to each other. In the reverse direction the backward traveling beam has a polarization orthogonal to the input polarizer and is therefore crossed with it, resulting in a rejected beam exiting the input polarizer.

1.2 Safe Use of Your Newport Faraday Isolator

The operational hazards presented to operating personnel by the use of your Newport Faraday Isolator are listed below. An explanation of how the Faraday Isolator is designed, together with procedures users can employ to eliminate or minimize these hazards are also listed.

1. Danger of sharp ferromagnetic objects being attracted to the residual permanent magnet fields outside of the isolator. This hazard is of most concern if such fields cause flying objects when being handled.

Your Newport Faraday Isolator requires strong internal magnetic fields to operate properly. Efforts have been made to minimize external fields from the device while still maintaining a relatively small and cost effective package. The external fields are designed to be well within federal safety guidelines which limit external fields from magnetic devices to be less than 2K Gauss at a radial distance of 5cm from the outside of the device. However, such fields can be sufficient to attract nearby objects such as knives and razor blades. Should attraction of such objects begin to occur there would be a strong attractive force directing these objects towards the interior of the magnet housing. This could be likely to result in injury (e.g. a cut or puncture wound) if such attraction occurred while the device was being handled – particularly if a body part of the operating personnel is near a beam Aperture (i.e. end) of the device.

To minimize the above risks remove all loose ferromagnetic objects from the path over which your Newport Faraday Isolator is to be moved prior to attempting to move it. **Do not** pick up the isolator by its ends (i.e. apertures) where the attractive magnetic fields are strongest. Always pick the isolator up along its sides.

2. Reflection of rejected beams from the input and output polarizer.

The polarizer covers have been positioned at the factory to block all beams rejected from the polarizers. In the event that your Faraday Isolator will be used with transmitted average powers in excess of 25W, or will block backward propagating light in excess of 0.5W average power, these polarizer covers must be rotated to allow rejected beams to exit (see Figure 1) onto user supplied beam dumps. **These rejected beams can represent a hazard to users and/or their colleagues. Care must be exercised to ensure that all rejected beams (both transmission and isolation directions) are accounted for and terminated into functional beam dumps.**

Wherever possible keep the strongest rejected beams in the horizontal plane of the table or otherwise safest direction (typically down into the table). **Always wear laser safety glasses/goggles consistent with all laser frequencies and power levels present.** See the following sections for further details.

3. Failure of operating personnel to observe standard laser safety by sighting down through the isolator when laser radiation is present.

The optical elements within the Newport Faraday Isolators can be transmissive throughout the visible and near infrared. Consequently it is never appropriate to view through the device in either the transmission or isolation direction when laser radiation is present – even with laser safety goggles.

Never sight through your Newport Faraday Isolator in either direction when there is any possibility of laser radiation being present.

4. Harm caused by external magnetic fields.

Your Newport Faraday Isolator has been designed to meet existing federal safety guidelines for external fields as noted previously. Such guidelines could change in the future as more information becomes known or reviewed regarding the interaction between magnetic fields and human health. Since there exist various claims regarding the potential harmful (and beneficial!) effects of magnetic fields on humans it is prudent to limit interaction with these fields as much as possible.

Personnel with any magnetically-sensitive implants such as pacemakers should consult their medical doctor regarding any potential complications which could arise from the isolator external magnetic fields.

5. Other non-health related hazards.

The Faraday Isolator external magnetic fields can draw ferromagnetic objects into the magnet housing that can damage the optical elements within the device. Keep a suitable area in all directions around the Faraday Isolator clear of any loose ferromagnetic objects. Ideally, use non-magnetic tools (such as stainless steel or titanium) and hardware to secure the Faraday Isolator. If only ferromagnetic tools are available use extreme care when using them around the Faraday Isolator. It is always helpful to bring such tools towards an aperture (or end) radially rather than along the optical beam path. Doing this ensures that the fields will tend to pull such objects into the magnet housing endplate rather than into the optical aperture. Where possible use two hands, one to hold the tool and the other to guide it to the desired destination.

Another concern regarding external magnetic fields is their effect on magnetically-sensitive devices. The external fields are strong enough to induce a pulse of current in electronic devices (such as digital watches) that can destroy them. The fields can also disrupt the operation of other mechanical devices with

ferromagnetic parts in them. Finally, the external fields can erase information from magnetic strips such as are found on credit and ID cards. Remove all magnetically-sensitive materials and devices such as watches, computer hard drives and magnetic strips from operators prior to working in the proximity of an isolator.

1.3 Operation

1.3.1 The Newport Broadband Faraday Isolator

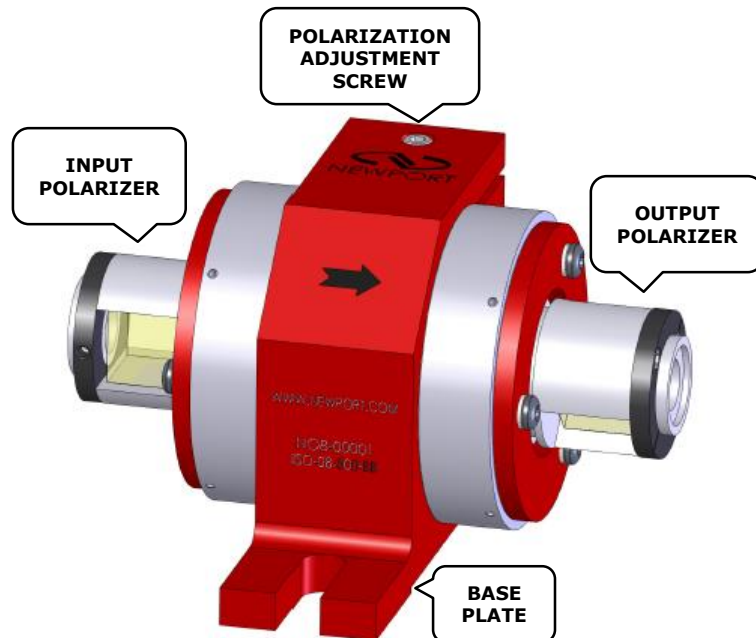


Figure 2: View of ISO-08-800-BB with polarizer covers removed

With the polarizer covers off, a polarizer can be seen at each end of the device. For both the 5mm and 8mm aperture devices, polarizing beam splitter cubes are used. The arrow on the side indicates the transmission direction. The input polarization shown is horizontal and the output is vertical. The central magnet housing together with the TGG and quartz crystals residing in its center form a broadband Faraday rotator. The input and output polarizers work in conjunction with the central Faraday rotator to form a Faraday Isolator as described previously in Section I. Figures 3 and 4 show a more detailed, diagram view of the device. Note that the input and output apertures are in-line and centered on the magnet body. This device may be adjusted readily for any input polarization. The polarization adjustment screw shown in Figure 2 may be loosened and the entire magnet housing along with input and output polarizers will rotate freely in the base clamp as a single assembly. Once the device is oriented for optimal

input polarization, the screw is tightened again. Further details for this procedure are provided in section 1.3.3.

The input view is observed in Figure 1. Note that the polarizer covers are in the open position here, allowing rejected beams to exit the device. The user may close these ports by simply rotating the polarizer cover, shown in Figure 1, until the port is in the closed position. This cover is held in place by an O-ring located between its inner surface and the polarizer mount.

Figure 3 shows a diagram with multiple views of the 5mm aperture Broadband Faraday Isolator.

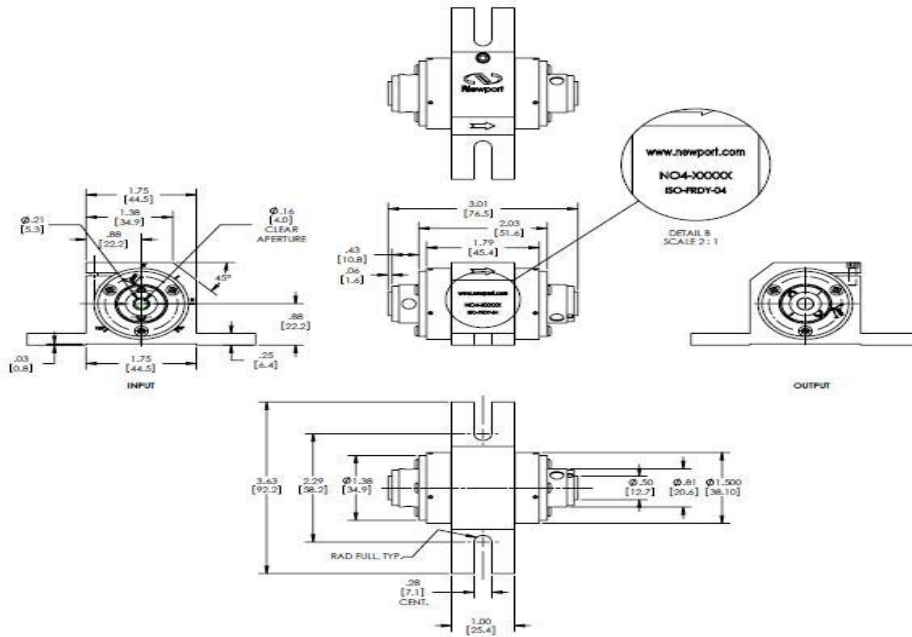


Figure 3: Diagram View of 5mm Aperture Broadband Isolator

Figure 4 shows a diagram with multiple views of the 8mm aperture Broadband Faraday Isolator.

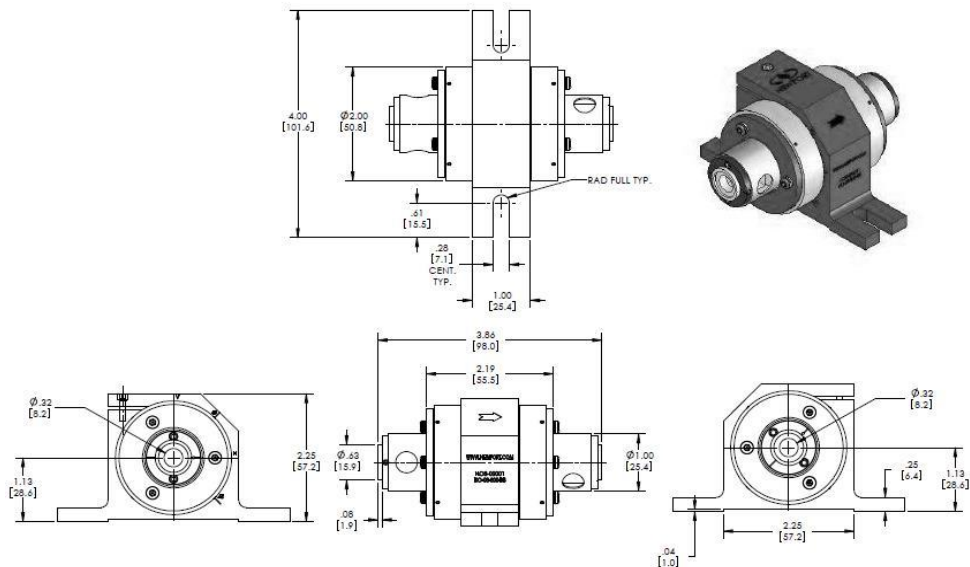


Figure 4: Diagram View of 8mm Aperture Broadband Isolator

1.3.2 Using your Faraday Isolator

Observe the guidelines for safe use of your Faraday Isolator found in Section 1.2 above when removing your isolator from its shipping container. Do not remove the protective dust-cover end caps from the polarizers until the device is in a clean, relatively dust-free environment. Save the protective end caps, packaging material and containers in the event that the device should ever need to be returned to Newport.

Verify that the Input and Output polarization states are consistent with the intended mode of operation. If not, re-adjust the isolator as required (see Section 1.3.3).

With the source laser off, or running at very low power (less than 100mW), position the Faraday Isolator such that the source laser beam can be directed through the Input Aperture.

Critical alignment of the Faraday Isolator should be done at low power (less than 100mW) in order to prevent optical damage to your isolator or laser source.

Use IR cards or viewers to ensure that the source laser beam is centered on the input and output apertures. The clear aperture of these devices is either 5mm or 8mm, centered on the circular cross section of the magnet body. There are different mounting options for establishing appropriate beam height. It is also preferable to use an IR viewer to ensure that weak reflections from AR coated optical surfaces in the Faraday Isolator are not being directed back into the source laser. The optical surfaces in the Faraday Isolator are angled slightly to reduce these reflections. Increasing the distance between the Faraday Isolator and the source laser can also help ensure that no reflections couple back into the source laser if necessary. Alternatively, if the beam used is smaller than the aperture of the device by a reasonable margin, the device may be slightly tilted.

At this point the Faraday Isolator should be secured to the work surface with two to four 1/4 - 20 or M6 screws – one for each slot in the baseplate flanges. Steel (ferromagnetic) ball drivers or other such wrenches will be attracted to the external magnetic field surrounding the device. If possible use anti-magnetic stainless steel or titanium tools. If ferromagnetic tools are used it is desirable to introduce them slowly toward the device from the sides along the direction of the baseplate flange slots.

If the Faraday Isolator will be used with average powers in excess of 25W transmitted or 0.5W rejected backward propagating radiation the Polarizer Covers will need to be removed so that the Escape Ports allow rejected polarization light to be safely dumped onto a beam dump. Failure to allow these rejected polarizations to escape can cause the device to heat up. Such heat can degrade the performance of the Faraday Isolator, or in severe cases, cause damage to optical components in the isolator. While working with low alignment level power and wearing safety glasses, remove the screw holding the dust cover in place for both the input and output. Any rejected polarized beams (in either the forward or backward propagating directions) can now exit the polarizer assembly. Use an IR viewer or IR card to locate these beams. Ensure that they are terminated on beam dumps consistent with the maximum amount of power that may be in such beams. If the Faraday Isolator is used in applications where strong reflections and/or optical gain elements (amplifiers) exist there may be very high power rejected

beams for backward propagating light at the input polarizer. If the average power levels used do not exceed 25W transmitted or 0.5W of backward propagating power then the Polarizer Covers may be kept in their factory positioned orientation – that is with all rejected beams blocked by the Polarizer Cover. However, if the Faraday Isolator is to be used with very high peak intensities it is prudent to allow rejected beams to escape on to external beam dumps to prevent any ablation damage to the nickel-plated Polarizer Covers. Follow the same procedure above as for high average powers in order to safely terminate all rejected beams.

1.3.3 Aligning your Faraday Isolator

For both the 5mm and 8mm aperture devices, the cylindrical magnet body is clamped into the base clamp structure. The polarization adjustment screw may be loosened and the cylindrical isolator structure rotates freely inside the clamp. The preferred method for alignment is to use an external polarizer mounted to a fine rotation stage with known directional axes and a waveplate to rotate the polarization before entering the test polarizer. By crossing this polarizer with the input polarizer of the device, a precise input polarization may be realized to a known reference orientation. (Note: the adjustment of the output polarizer must be 90 degrees from the input polarizer.) For more details, please contact Newport’s technical sales department at tech@newport.com.

1.4 Specifications

Model	ISO-05-800-BB	ISO-08-800-BB
Polarization Direction of Input Beam	Horizontal Polarization	Horizontal Polarization
Clear Aperture (mm)	5	8
Center Wavelength (nm)	800	800
Spectral Range (nm)	720-950	720-950
Isolation (dB) @ 22°C	>30	>30
Transmission (%) @ 22°C	>92	>92
Pulse Damage Threshold (J/cm ²) For a 10ns pulse	1	1
Polarizer Type	PBS Cube	PBS Cube
Rotating Medium	Terbium Gallium Granite	Terbium Gallium Granite

Notes:

1. Polarization Orientation: Please note that the Broadband Isolators utilize a 45° crystal quartz rotator. At the center wavelength there will be 90° polarization rotation in the forward direction. In the reverse direction, there will be 0° polarization rotation across the entire spectral bandwidth of the device. For more information on how our Broadband Isolators rotate the plane of polarized light, please view our Application Notes available through the Support tab.

2. Dispersion: Some pulse broadening does occur when using Newport's Broadband Isolator. Use the following Sellmeier Equation:

$$n^2 - 1 = \frac{E_d E_o}{E_o^2 - (\hbar c / \lambda)^2}$$

where: $E_o = 9.223\text{eV}$ and $E_d = 25.208\text{eV}$

3. Operating and Storage Temperature: Performance of Newport's Broadband Isolators is related to operating temperature.

2

Factory Service Information

2.1 Service Form



Newport Corporation
U.S.A. Office: 800-222-6440
FAX: 949/253-1479

Name _____ Return Material Authorization # _____
(Please obtain RMA# prior to return of item)

Company _____

Address _____ Date _____

Country _____ Phone Number _____

P.O. Number _____ FAX Number _____

Item(s) Being Returned:

Model # _____ Serial # _____

Description

Reason for return of goods (please list any specific problems)

**Newport Corporation
Worldwide Headquarters**

1791 Deere Avenue
Irvine, CA 92606

(In U.S.): 800-222-6440
Tel: 949-863-3144
Fax: 949-253-1680

Internet: sales@newport.com



Newport
Experience | Solutions

Visit Newport Online at: www.newport.com



Newport Corporation, Irvine, California, has
been certified compliant with ISO 9001 by
the British Standards Institution.

Printed in the U.S.