


# Compact Isolators

## Instruction Manual

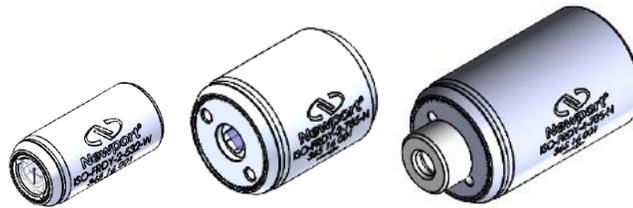

Optical Isolators  
405 nm to 1064 nm


Thank you for purchasing your **Compact** Optical Isolator from Newport. This user guide will help answer questions you may have regarding the safe use and optimal operation of your Optical Isolator.

## TABLE OF CONTENTS


I. Compact Optical Isolator Overview .....	2
II. Safe use of your Compact Optical Isolator .....	3
III. The Compact Optical Isolator .....	4
IV. Using your Compact Optical Isolator .....	5
V. Tuning your Compact Optical Isolator .....	6
VI. Warranty Statement and Repair .....	6

### I. Compact Optical Isolator Overview

Your **Compact** Optical Isolator is essentially a unidirectional light valve. It is used to protect a laser source from destabilizing feedback or actual damage from back-reflected light. Figure 1 and 2 below identifies the main elements of your Optical Isolator.

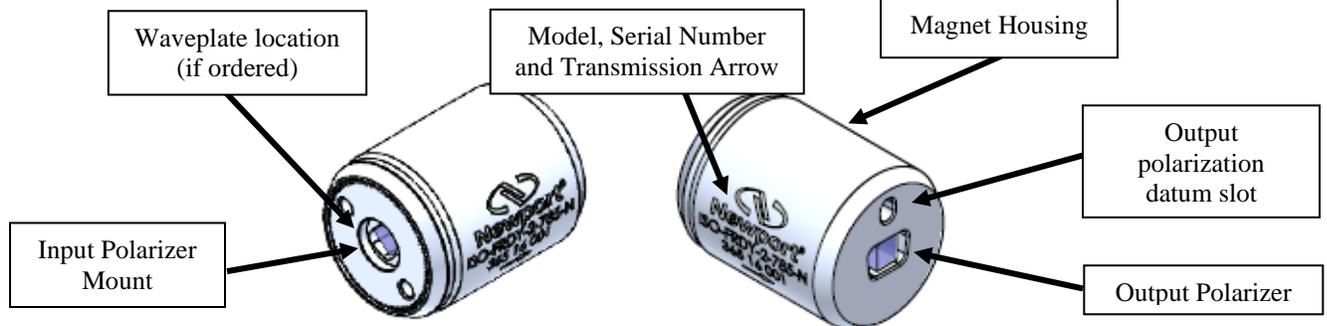

**Figure 1: 2 mm Compact Optical Isolator**

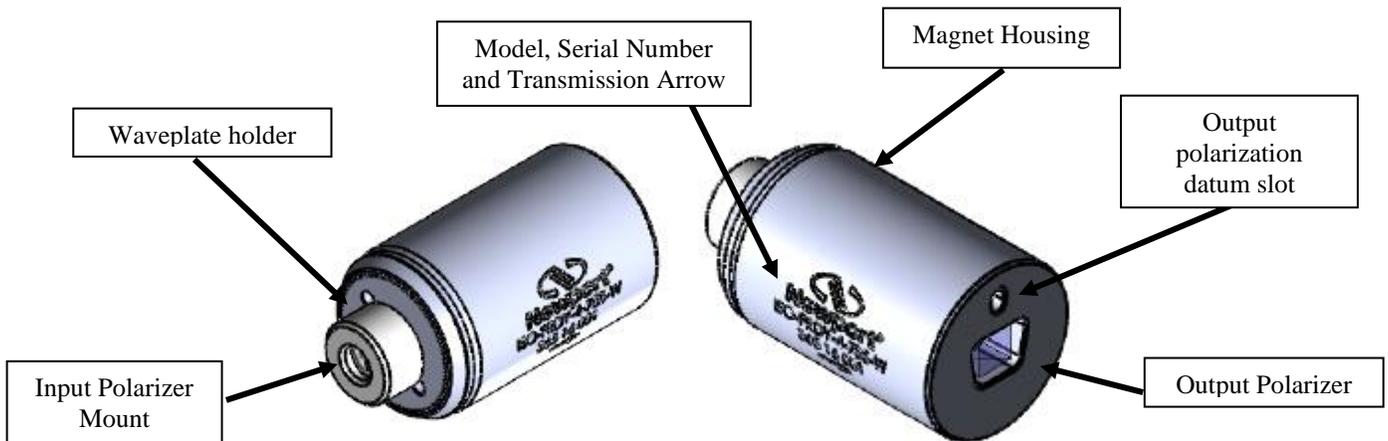

**Figure 2: 4 mm Compact Optical Isolator**

The Optical Isolator is a cylindrically-shaped magneto-optic device. Strong Neodymium Iron Boron permanent magnets are used to generate high (>10,000 Gauss) axially-

oriented fields within the magnet housing. The strong longitudinal field causes 45 degrees of non-reciprocal polarization rotation for propagating light via the Faraday Effect in the terbium gallium garnet (TGG) crystal located within the magnet housing. In operation, the magnet housing is sandwiched between input and output polarizers that have their transmission axis oriented 45 degrees relative to each other to account for the 45 degrees of Faraday rotation in the TGG crystal in the forward (transmission) direction. In the reverse (isolation) direction, the non-reciprocal Faraday rotation and the 45 degree polarizer transmission axis angle add so that the polarization transmitted by the output polarizer is rejected at the input polarizer.

Your **Compact** Optical Isolator is labeled with a serial number on the magnet body of the device.

## II. Safe use of your Compact Optical Isolator

The operational hazards presented to operating personnel by the use of your **Compact** Optical Isolator are listed below. An explanation of how the Optical Isolator is designed together with procedures users can employ to eliminate or minimize these hazards is presented.

1. Danger of sharp ferromagnetic objects being attracted to the residual permanent magnet fields outside of the Faraday Rotator. This hazard is of most concern if such fields cause flying objects when being handled.

Your Isolator requires strong internal magnetic fields to operate properly. Efforts have been made to minimize external fields from the device while still maintaining a relatively small and cost-effective package. The external fields are designed to be well within federal safety guidelines which limit external fields from magnetic devices to be less than 2 K Gauss at a radial distance of 5cm from the outside of the device. However, such fields can be sufficient to attract nearby objects such as knives and razor blades. Should attraction of such objects begin to occur, there would be a strong attractive force directing these objects towards the interior of the magnet housing. This could be particularly likely to result in injury (e.g. a cut or puncture wound) if such attraction occurred while the device was being handled – particularly if a body part of the operating personnel is near a beam aperture (i.e. end) of the device.

To minimize the above risks, remove all loose ferromagnetic objects from the path over which your Compact Optical Isolator is to be moved prior to attempting to move it. **Do not** pick up the isolator by its ends (i.e. apertures) where the attractive magnetic fields are strongest. Always pick the isolator up along its sides.

2. Never attempt to disassemble the magnetic housing of your Faraday Isolator. Serious injury could result.
3. Failure of operating personnel to observe standard laser safety by sighting down through the Faraday Rotator when laser radiation is present.

The optical elements within the Compact Optical Isolators can be transmissive throughout the visible and near infrared. Consequently, it is never appropriate to view through the device in either the transmission or isolation direction when laser radiation is present – even with laser safety goggles.

***Never sight through your Optical Isolator in either direction when there is any possibility of laser radiation being present.***

4. Harm caused by external magnetic fields.

Your **Compact** Optical Isolator has been designed to meet existing federal safety guidelines for external fields as noted previously. Such guidelines could change in the future as more information becomes known or reviewed regarding the interaction between magnetic fields and human health. Since various claims exist regarding the potential harmful (and beneficial!) effects of magnetic fields on humans it is prudent to limit interaction with these fields as much as possible.

Personnel with any magnetically-sensitive implants such as pacemakers should present a copy of this report and consult their medical doctor regarding any potential complications which could arise from the isolator external magnetic fields.

5. Other non-health related hazards.


The Optical Isolator external magnetic fields can draw ferromagnetic objects into the magnet housing which can damage the optical elements within the device. Keep a suitable area from the Optical Isolator in all directions clear of any loose ferromagnetic objects. Ideally, use non-magnetic tools (such as stainless steel or titanium) and hardware to secure the Optical Isolator. If only ferromagnetic tools are available use extreme care when using them around the Optical Isolator. It is always helpful to bring such tools towards an aperture (or end) radially rather than along the optical beam path. Doing this ensures that the fields will tend to pull such objects into the magnet housing endplate rather than into the optical aperture. Where possible use two hands, one to hold the tool and the other to guide it to the desired destination.

Another concern regarding external magnetic fields is their effect on magnetically-sensitive devices. The external fields are strong enough to induce a pulse of current in electronic devices (such as digital watches) that can destroy them. The fields can also disrupt the operation of other mechanical devices with ferromagnetic parts in them. Finally, the external fields can erase information from magnetic strips such as those found on credit and ID cards. Remove all magnetically-sensitive materials and devices such as watches, computer hard drives and magnetic strips from operators prior to working in the proximity of an isolator.


### III. The Compact Optical Isolator

A polarizing beamsplitter cube (PBSC) can be seen inside the input polarization mount. The inscribed arrow on the magnet body shows the transmission direction. The output

PBSC is seen to be oriented with its transmission axis rotated 45 degrees relative to the input PBSC. The input polarization shown is horizontal. The central magnet housing together with the TGG crystal residing in its center forms a Faraday Rotator. The Faraday Rotator rotates the input horizontal transmission axis by 45 degrees so that transmitted light has a polarization aligned with the output transmission axis. The input and output PBSCs work in conjunction with the central Faraday Rotator to form an Optical Isolator as described previously in Section I. Though the overall size of the device varies depending on the wavelength dependent model, the size and operation of the PBSC mounts are identical. Figures 3 below shows output views of a **Compact Optical Isolator**. Figure 4 shows polarization orientation reference.

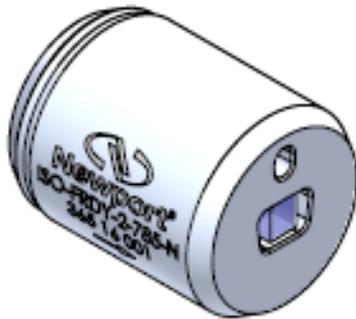

**Figure 3: Horizontal Input Polarization**

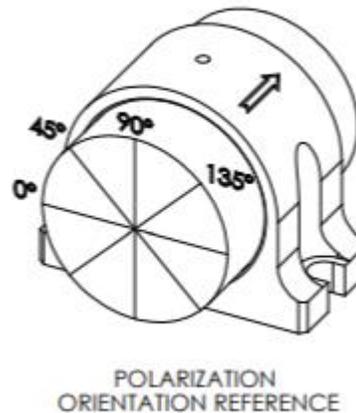

**Figure 4: Polarization Reference**

The polarization (vertical or horizontal) is controlled by the slot on the output of the device. When the locating slot is vertical as shown in Figure 3, the output polarization is 90 degrees. This can be changed by simply rotating the device such that the slot is angled to that of the desired output polarization. The cubes cannot be adjusted, they are set at the factory and the rotation is optimized to 45 degrees. These principals are valid with all sizes of **Compact Optical Isolators**. For all 4 mm devices, the waveplate holder will be included on your device whether a waveplate is ordered or not. For the 2 mm devices, there will be no waveplate holder present, whether or not a waveplate is ordered.

#### IV. Using your Compact Optical Isolator

Observe the guidelines for safe use of your Optical Isolator found in Section II above when removing your isolator from its shipping container. Do not remove from shipping container until the device is in a clean, relatively dust-free environment. Save the packaging container in the event that the device should ever need to be returned to Newport.

With the source laser off, or running at very low power (less than 250 mW), position the Optical Isolator such that the source laser beam can be directed through the input

aperture. Critical alignment of the Optical Isolator should be done at low power (less than 250 mW) in order to prevent optical damage to your isolator or laser source.

Use IR cards, IR viewers, or heat-sensitive film to ensure that the source laser beam is centered on the input and output apertures. It is also preferable to use an IR viewer to ensure that weak reflections from AR coated optical surfaces in the Optical Isolator are not being directed back into the source laser. The optical surfaces in the Optical Isolator are angled slightly to reduce these reflections. Increasing the distance between the Optical Isolator and the source laser can also help ensure that no reflections couple back into the source laser if necessary.

At this point, the Optical Isolator should be secured into the final location.

If the Optical Isolator will be used with average powers in excess of 12 W transmitted or 0.5 W rejected backward propagating radiation please consult damage threshold specifications for operating range. Since there is no way for the rejected polarizations to escape, at higher power it can cause the device to heat up. Such heat can degrade the performance of the Optical Isolator, or in severe cases, cause damage to optical components in the isolator.

## V. Tuning your Compact Optical Isolator

### A. *Adjusting Input Polarization*

The transmission axis of the output polarizer is 45 degrees from the slot position on the magnet body on the output side of the device. If the linear polarization of the laser source is geometrically known, aligning the input polarization of the Optical Isolator to that of the laser source is straightforward. Simply rotate the magnet housing until the input polarizer transmission axis is aligned to that of your laser source. Wear laser safety glasses or goggles. Protect your colleagues by ensuring that they are wearing safety glasses or goggles as appropriate.

### B. *Fine Wavelength Adjustment*

The wavelength of the device is determined at time of order and optimized to that during the device assembly. The device cannot be tuned by the customer. Please consult the data sheet to usable wave ranges.

### C. *Waveplate Option and Adjustment*

The waveplate is installed at the factory and cannot be adjusted. Should you need to adjust the waveplate please contact.

## VI. Warranty Statement and Repair

## Warranty

Newport Corporation warrants that this product will be free from defects in material and workmanship and will comply with Newport's published specifications at the time of sale for a period of one year from date of shipment. If found to be defective during the warranty period, the product will either be repaired or replaced at Newport's option.

To exercise this warranty, write or call your local Newport office or representative. You will be given prompt assistance and return instructions. Send the product, freight prepaid, to the indicated service facility. Repairs will be made, and the instrument returned freight prepaid. Repaired products are warranted for the remainder of the original warranty period or 90 days, whichever first occurs.

### Limitation of Warranty

The above warranties do not apply to products which have been repaired or modified without Newport's written approval, or products subjected to unusual physical, thermal or electrical stress, improper installation, misuse, abuse, accident or negligence in use, storage, transportation or handling. This warranty also does not apply to fuses, batteries, or damage from battery leakage.

THIS WARRANTY IS IN LIEU OF ALL OTHER WARRANTIES, EXPRESSED OR IMPLIED, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR USE. NEWPORT CORPORATION SHALL NOT BE LIABLE FOR ANY INDIRECT, SPECIAL, OR CONSEQUENTIAL DAMAGES RESULTING FROM THE PURCHASE OR USE OF ITS PRODUCTS.

## Confidentiality & Proprietary Rights

### Reservation of Title:

The Newport programs and all materials furnished or produced in connection with them ("Related Materials") contain trade secrets of Newport and are for use only in the manner expressly permitted. Newport claims and reserves all rights and benefits afforded under law in the Programs provided by Newport Corporation.

Newport shall retain full ownership of Intellectual Property Rights in and to all development, process, align or assembly technologies developed and other derivative work that may be developed by Newport. Customer shall not challenge or cause any third party to challenge the rights of Newport.

### Preservation of Secrecy and Confidentiality and Restrictions to Access:

Customer shall protect the Newport Programs and Related Materials as trade secrets of Newport and shall devote its best efforts to ensure that all its personnel protect the Newport Programs as trade secrets of Newport Corporation. Customer shall not at any time disclose Newport's trade secrets to any other person, firm, organization, or employee that does not need (consistent with Customer's right of use hereunder) to obtain access to the Newport Programs and Related Materials. These restrictions shall

not apply to information (1) generally known to the public or obtainable from public sources; (2) readily apparent from the keyboard operations, visual display, or output reports of the Programs; (3) previously in the possession of Customer or subsequently developed or acquired without reliance on the Newport Programs; or (4) approved by Newport for release without restriction.

First printing 2019

© 2019 by Newport Corporation, Irvine, CA. All rights reserved.

No part of this manual may be reproduced or copied without the prior written approval of Newport Corporation.

This manual has been provided for information only and product specifications are subject to change without notice.

Any change will be reflected in future printings.

Newport Corporation 1791 Deere Avenue Irvine, CA, 92606 USA

## Technical Support Contacts

### North America

#### Newport Corporation

1791 Deere Ave, Irvine, CA 92606

Telephone: (877) 835-9620

Telephone: (949) 863-3144

### Europe

Newport/MICRO-CONTROLE S.A.

Zone Industrielle

45340 Beaune la Rolande,

FRANCE

Telephone: (33) 02 38 40 51 56

### Asia

#### Newport Opto-Electronics Technologies

中国 上海市 爱都路 253号 第3号楼 3层  
C部位, 邮编 200131

253 Aidu Road, Bld #3, Flr 3, Sec C,  
Shanghai 200131, China

Telephone: +86-21-5046 2300

Fax: +86-21-5046 2323

### Newport Corporation Calling Procedure

If there are any defects in material or workmanship or a failure to meet specifications, promptly notify Newport's Returns Department by calling 1-800-222-6440 or by visiting our website at [www.newport.com/returns](http://www.newport.com/returns) within the warranty period to obtain a **Return Material Authorization Number (RMA#)**. Return the product to Newport Corporation, freight prepaid, clearly marked with the RMA# and we will either repair or replace it at our discretion. Newport is not responsible for damage occurring in transit and is not obligated to accept products returned without an RMA#.

E-mail: [rma.service@newport.com](mailto:rma.service@newport.com)

When calling Newport Corporation, please provide the customer care representative with the following information:

- Your Contact Information
- Serial number or original order number
- Description of problem (i.e., hardware or software)

To help our Technical Support Representatives diagnose your problem, please note the following conditions:

- Is the system used for manufacturing or research and development?
- What was the state of the system right before the problem?
- Have you seen this problem before? If so, how often?
- Can the system continue to operate with this problem? Or is the system non-operational?
- Can you identify anything that was different before this problem occurred?

### Notes:

**Operating Temperature:** Performance of Newport's Faraday Rotators/Isolators is related to

operating temperature. For information on the effect of operating temperature on Newport's Faraday rotators/isolators, please review our application note, *Thermal Lensing Analysis of TGG*.

**For higher incident powers or aperture sizes other than those specified, or other custom requirements, please consult Newport.** You may contact Newport's sales department at [sales@newport.com](mailto:sales@newport.com).